

User Instruction Manual


Hoyer® Slings

To avoid injury, read user manual prior to use.


MODELS

Available in solid poly (P)
and poly mesh fabrics (D):


9 Commode seat,
back and set of
two chains

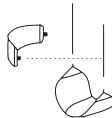

105 One piece XL solid
and mesh commode
seat and back with
set of two chains.

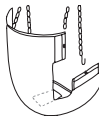


110 Comfort seat with
wire insert head
support and set of
two web straps

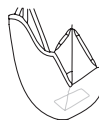

111 Comfort commode
seat with wire
insert head support
and set of two
web straps




112 One piece seat and
back with set of two
chains

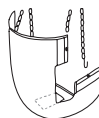


113 One piece commode
seat and back with
set of two chains

114 Comfort seat, less
head support and
set of two web
straps

115 Comfort commode
seat, less head
support and set of
two web straps

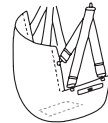


117 Comfort commode
seat with modified
head support and
set of two chains

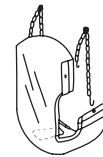

118 Comfort seat with
modified head
support and set of
two chains

Available in solid poly only (P):

127 Split commode seat
sling with wire insert
head support and set
of two web straps

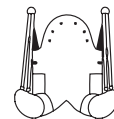

DIRECTIONS:

1. Arrange all items for a smooth transfer - lifter, sling, wheelchair.
2. Patient should be in center of twin size bed or one side of double size bed.
3. **Models 110, 111, & 127 only** - Insert head wire into pockets of the back part of the comfort seat). Slip head piece of sling over head wire. Close snaps.
4. Roll patient on side away from attendant. (raise side rail, if equipped, on side opposite attendant).
5. **Model 127 only** - Place sling folded half-way under patient so cut-out is just below tail-bone.
6. **Model 9 only** - Place wider piece (seat) under patient's thighs so lower edge of seat is up to knees. Place narrower piece (back) just above small of back. If patient is in hospital bed, position seat sling, then elevate head of bed to facilitate placing back piece.
7. **All other models** - Place sling folded half-way under patient so lower edge of seat is slightly below the knees.
8. Roll patient towards attendant. Pull sling through (like positioning draw sheet).
9. **Model 127 only** - If patient is in chair, you may place sling under patient without lifting patient by following these few simple steps:
 - a. Have Leg flaps open.
 - b. Have patient lean forward slightly and slide open portion of sling down to seat and lay leg flaps alongside patient's legs.
 - c. Grasp leg flap on one side of patient and, while holding against knee, pull leg flap forward; repeat with other side until sling is in correct position.
 - d. Bring flap under one thigh and insert ring (A) into snap (A). Repeat for other leg, inserting ring (B) into snap (B), thus enclosing a thigh in each leg flap.
 - e. If you wish to use the sling without containing the thighs, follow instructions above but bring ring (A) to snap (B) and ring (B) to snap (A), criss-crossing flaps under the leg.
10. It will help to raise the head of the bed after sling in positioned if patient is in a hospital bed.
11. **STEPS 11-15 FOR MODEL 9 ONLY** - Attach the S-hooks of the chain to the loops of the seat hangers. Hooks should be inserted AWAY from the patient to the outside of the sling.
12. Attach the ends of the chains to the swivel bar hooks. You may hook any link of chain into the swivel bar hook to adjust the heights of sling from floor. Make sure links are equal on each side.
13. Attach the S-hooks of the back in links 1,2 or 3, as required. Count links to be sure there are the same number on each side. Check to see that the S-hooks are hooked all the way into the chains. **DO NOT PLACE S-HOOKS OF THE BACK INTO EXCESS PORTION OF CHAIN BUT INTO LINKS BETWEEN SEAT AND SWIVEL BAR!**
14. Patient's arms should be outside chains if possible. They may hold on the chains if desired.
15. Check chains and S-hooks to see that they are properly positioned.
16. **STEPS 16-19 OTHER SLINGS WITH CHAIN ONLY.** The 133-S-C chain has the 9th link painted red. Attach the S-Hook closest to this link into the hole of the back part. Attach the S-Hook of the opposite end into the hole in the seat part of the sling (one chain attached per side). Make sure the S-hooks are hooked all the way into the holes. Hooks should be inserted AWAY from the patient and to the outside of sling.
17. Move the lifter so open of U-base or end of C-base is under side of bed.

18. Hook the 9th link of the chain (painted red) into the end of the swivel bar. Check to see that links are hooked all the way into the swivel bar. Posture of patient can be changed by varying the link hooked into the swivel bar. It is important that the section of chain from the swivel bar to the patients back be SHORTER than the section between swivel bar and thigh of patient.
19. Patient's arm may be outside of sling, if desired.
20. **STEPS 20-24 FOR SLINGS WITH WEB STRAPS ONLY.** S-Hook of the shortest straps of the web straps are hooked into the holes of the back. Make sure S-Hooks are hooked all the way into the holes. Hooks should be inserted AWAY from the patient and to the outside of the sling. Hole to use varies for each individual. Practice to see which is more comfortable.
21. Hook the center strap into the top hole of the seat and the last strap into the hole of the seat. Make sure S-hooks are hooked all the way into the holes. Hooks should be inserted AWAY from the patient to the outside of the sling.
22. Move the lifter so open end of U-base or end of C-base is under side of bed.
23. Hook the D-rings of the webs straps into the swivel bar hooks.
24. Check to see that lower edge of seat is still close to knees and S-hooks are still properly positioned.
25. Pump hydraulic jack handle. Partially lift the patient; check S-hooks and swivel bar for proper hook up. Continue pumping. Steering handle may be held to steady the pumping.
26. Raise patient until buttocks are just clear of the bed. Swing patient's feet off bed and turn patient to face mast. Grasp steering handles and move patient away from bed. Lower patient as soon as they are clear of bed to keep center of gravity low. Position patient over chair or commode. You may have to raise them from the lower position. U-base or C-base fits around chair. Wheelchair may be brought and positioned under patient. Lock wheelchair locks prior to lowering patient into chair.
- 27a. If the jack has a pin release, press handle SLOWLY towards pump body. Lower patient slowly. Guide their decent. Push gently on their knees as they are being lowered so correct sitting position will result. After reaching seat, press handle against pin and press down on boom.
- 27b. If the jack has to turn now release, lower patient by turning the release knob gently to the left. Guide their decent. Push gently on their knees as they are being lowered so correct sitting position will result. After reaching seat, open release knob a turn and press down on boom.
28. Detach web straps/chains from swivel bar and move lifter away. Then remove web straps/chains from sling. Patient may remain on sling or sling may be removed in reverse of above.

TO RETURN PATIENT TO BED:

1. Same procedure in reverse. Check to see that S-hooks are properly positioned in sling. Make sure S-hooks are not caught on chair arms as patient is being lifted.
2. Center patient over bed and lower gently.

Continued on rear page >

DIRECTIONS:

1. Keep web straps or chains taut to hold patient in a comfortable sitting position.
2. Adjust clothing before moving lifter to straddle commode.

USE CARE, DISCRETION AND COMMON SENSE TO DETERMINE IF A SEVERELY DISABLED OR HANDICAPPED PERSON CAN BE LIFTED WITH A HOYER LIFTER.

PLEASE NOTE: Web straps (W-C) are the only set to be used with the 110, 111, 114, 115 and 127 slings. 133-S-C or 132-S-C chains are the only set to be used with 112, 113, 117, 118 and 123 slings. 134-S-C chains are the only chain set to be used with the 9 sling.

WASHING INSTRUCTIONS:

Machine wash warm 104°F (40°C) or cold. Do not exceed 185°F (85°C) water temperature. Air dry or very low dryer heat. **DO NOT USE BLEACH.** Do not wash with other colours. Remove bars before washing. (The bars and seat hangers of the two-piece slings are not removable. It is recommended that these be hand washed).

 WARNING

Always check slings before use. Bleached, torn, cut frayed or broken slings are unsafe and could result in injury. Destroy and discard. Do not alter slings.

 WARNING

Do not wash slings with bleach. Bleach will damage material. Slings can suffer damage during washing and drying.

Questions regarding the proper care and use of your sling may be directed to Customer Service at 1-800-826-0270.