




ProCair Plus
and Optional P.R.O. Technology® ES2 Control Unit
Convertible Alternating Pressure Therapy System

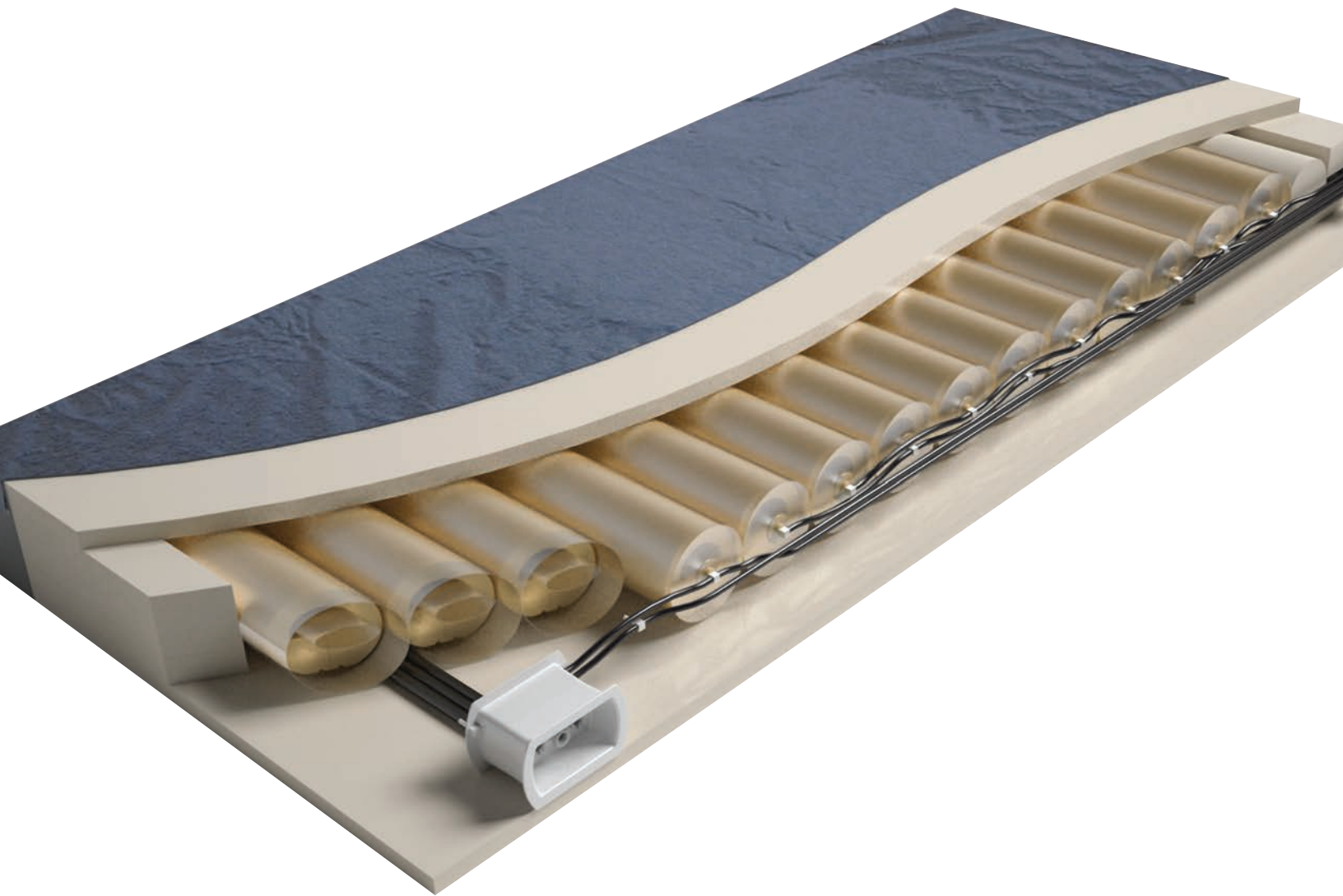

transforming wound management

ProCair Plus and Optional Control Unit

Non-Powered, Convertible Alternating Pressure Therapy System



Features and Benefits

Mattress Therapy Features

- Three-zoned system (shoulder, torso and foot) inner air core automatically adjusts to each patient's body profile and weight
- Fifteen cells with 14 providing alternating pressure therapy with optional control unit
- P.R.O. Technology valve maintains optimal air pressures in each therapy zone for superior pressure redistribution regardless of patient weight and profile
- One-inch viscoelastic topper enhances pressure redistribution and patient comfort
- Surrounding firm perimeter provides stability during patient care and contributes to patient safety by eliminating gaps between frames and siderails
- Therapy pad with these unique qualities:
 - Fluid proof and stain resistant
 - Low friction and low shear surface
 - Zippered construction
- Available in multiple widths and lengths
- Meets flammability standards TB 117, TB 129, 16 CFR 1633 and BFD IX-11



Optional P.R.O. Technology® ES2 Control Unit

- Adds your ProCair Plus to alternating pressure therapy to enhance patient comfort and wound healing
- Three alternating pressure zones, including the delicate heel section for optimal pressure redistribution
- Three modes of operation – *Autofirm*, *Therapy* and *Alternating*
 - *Autofirm* assists patients and caregivers during patient transfer and treatment
 - *Therapy* control settings provides powered immersion and flotation for excellent pressure redistribution.
 - *Alternating* cycle time options of 5, 10 and 15 minutes
- Closed loop control system eliminates concerns of changes in mattress interface pressure due to patient repositioning and ambient temperature changes

Clinical Indications

Pressure Relief: Stage I – II Pressure Ulcers¹; Rehabilitation; Burns; Neurology; Dermatology; Amputations

Pain Management: AIDS; Arthritis; Oncology

¹Stage II in non-powered mode and up to uncomplicated Stage III or IV using alternating pressure mode with optional control unit based on individual patient assessment.

Technical Specifications

Maximum Weight Capacity*: 500 lbs

Control Unit

Weight: 7 lbs.

Dimensions: 7.5" (W) x 12.25"(H) x 5.5"(D)

Mattress Dimensions

36"(W) x 80"(L) x 7"(D)

36"(W) x 84"(L) x 7"(D)

42"(W) x 80"(L) x 7"(D)

Agency Approvals (ES Control Unit)

UL Classified: UL 60601-1 CAN/CSA C22.2 No.60601-1

Electrical Specifications (Control Unit)

120V AC, 60 Hz, 0.6A

Environmental Conditions

Operating Conditions

Ambient Temperature: +10°C to +40°C

Relative Humidity: 30% to 75% Non-Condensing

Storage and Shipping Conditions

Ambient Temperature: -10°C to +40°C

Relative Humidity: 10% to 100%

*Mattress weight capacity only; total weight must not exceed bed frame manufacturer's specified load capacity.

Joerns reserves the right to change specifications without notice.



Post Acute, Acute, HomeCare

2430 Whitehall Park Drive
Charlotte, NC 28273
(P) 800.826.0270
(F) 800.457.8827

VA/Government

19748 Dearborn Street
Chatsworth, CA 91311
(P) 800.966.6662
(F) 800.232.9796

Canadian Office

1000 Clarke Road, Suite 6
London, ON Canada N5V 3A9
(P) 866.546.1151
(F) 519.451.8662