




ReplaMatt Basic

Therapeutic Mattress Replacement

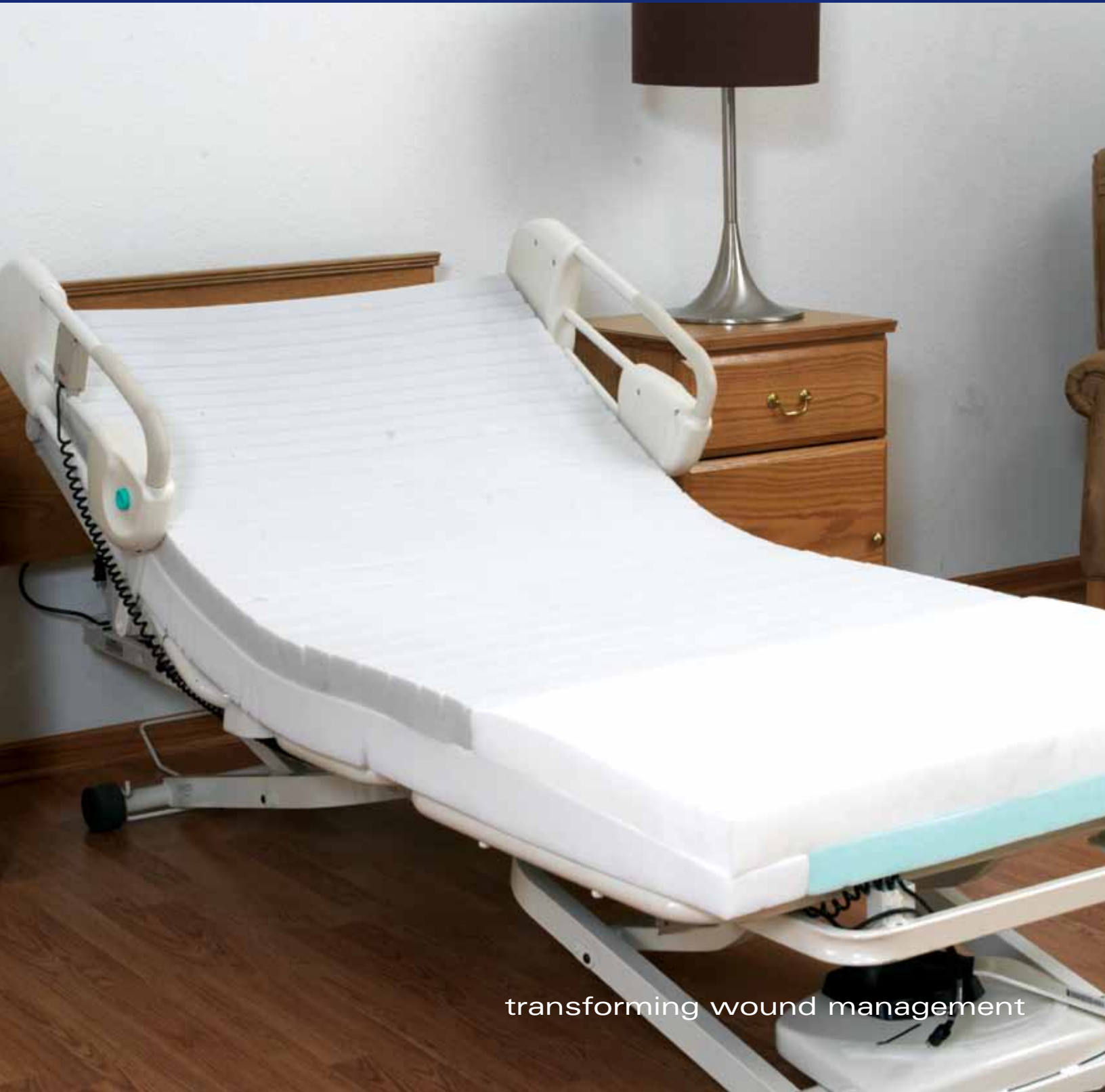

transforming wound management



Features and Benefits

Mattress Therapy Features

- Firm, yet highly resilient foam, gently cradles head and torso for maximum patient comfort
- Integrated viscoelastic foam heel slope provides a gentle tapering of the mattress that redistributes pressure along the lower legs while reducing pressure in the delicate heel section
- Strategically located articulation cuts in mattress base allow the mattress to move in harmony with the bed frame reducing friction and shear on patient skin and increasing product life cycle
- Firm perimeter and Tru-fit sizing help meet FDA/HBSW guidelines
- 70-denier therapy cover with these unique qualities:
 - Fluid proof and stain resistant
 - Infection control treatment
 - Low friction and low shear surface
- Non-skid bottom prevents mattress sliding on the deck to reduce the risk of falls, is tear resistant
- Latex free mattress construction

Clinical Indications

Pressure Relief: Recommended for prevention and treatment of pressure ulcers up to Stage I.

Technical Specifications

Maximum Weight Capacity*: 500 lbs

Mattress

Dimensions: 35"(W) x 80"(L) x 7"(D)

36"(W) x 80"(L) x 7"(D)

36"(W) x 84"(L) x 7"(D)

Warranty

Mattress: Five-years

Cover: One-year

Flammability Standards

Fire barrier meets fire standards TB 117, TB 129,
16 CFR 1633 and BFD IX-11

ILD and Density

Top Layer – 1.7 lb. / 15 ILD

Top Layer Heel – 3.0 lb. / 15 ILD

Perimeter – 1.7 lb. / 55 ILD

Base Layer – 1.7 lb. / 30 ILD

**Mattress weight capacity only; total weight must not exceed bed frame manufacturer's specified load capacity.*
Joerns reserves the right to change specifications without notice.