




# ReplaMatt Deluxe

Therapeutic Mattress Replacement

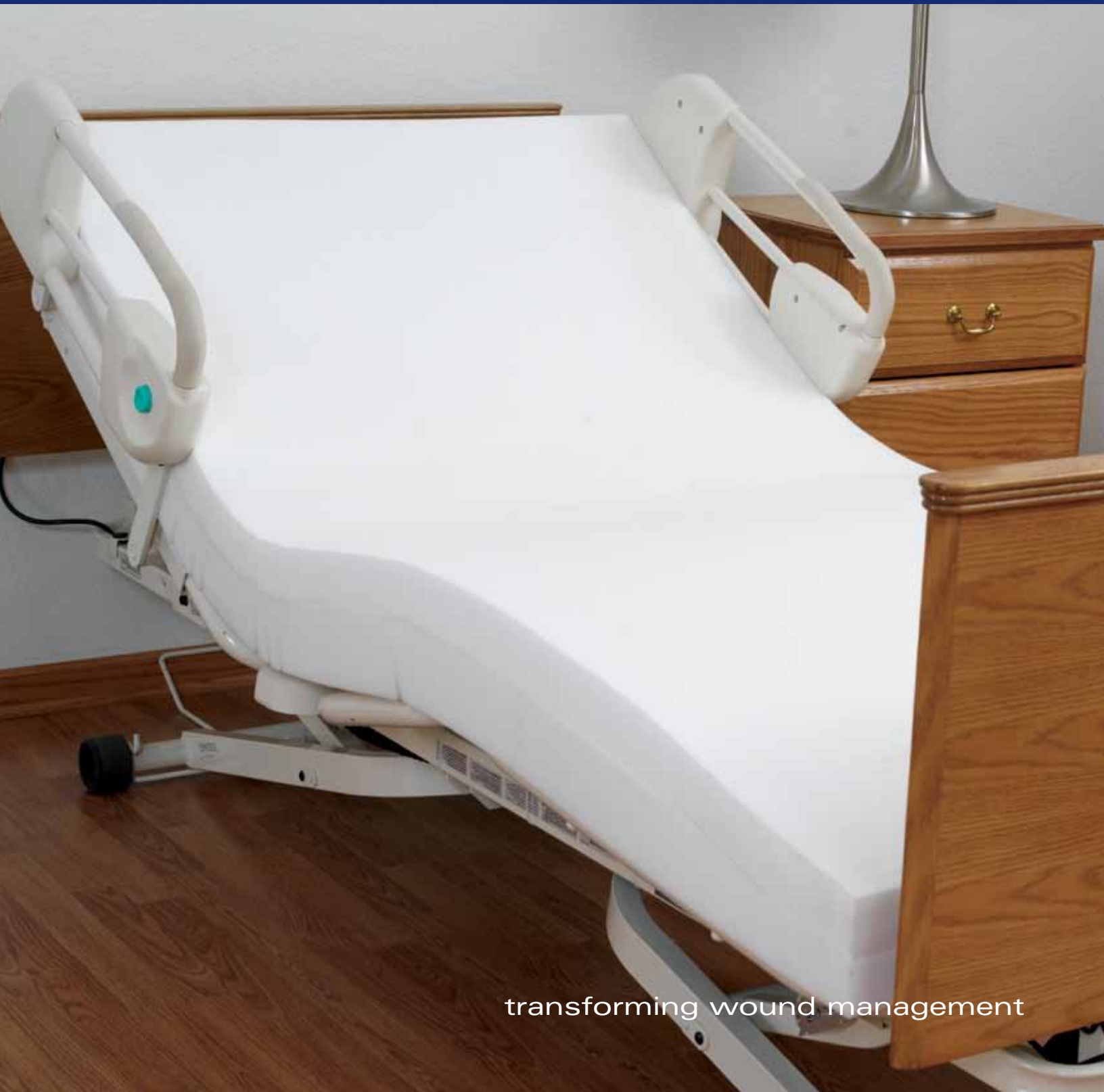

transforming wound management



## Features and Benefits

### Mattress Therapy Features

- Two-ply construction incorporates cross-cut two-inch viscoelastic foam topper to gently cradle head and torso for maximum patient comfort
- Integrated heel slope provides a gentle tapering of the mattress that redistributes pressure along the lower legs while reducing pressure in the delicate heel section
- Strategically located articulation cuts in mattress base allow the mattress to move in harmony with the bed frame reducing friction and shear on patient skin and increasing product life cycle
- Firm perimeter and Tru-fit sizing help meet FDA/HBSW guidelines
- 70-denier therapy cover with these unique qualities:
  - Fluid proof and stain resistant
  - Infection control treatment
  - Low friction and low shear surface
- Tear resistant, non-skid bottom prevents mattress sliding on the deck to reduce the risk of falls
- Latex free mattress construction

## Clinical Indications

**Pressure Relief:** Recommended for prevention and treatment of up to Stage II pressure ulcers.

## Technical Specifications

**Maximum Weight Capacity\*:** 500 lbs

### Mattress

Dimensions: 35"(W) x 80"(L) x 7"(D)

36"(W) x 80"(L) x 7"(D)

36"(W) x 84"(L) x 7"(D)

### Warranty

Mattress: Five-years

Cover: One-year

### Flammability Standards

Fire barrier meets fire standards TB 117, TB 129,  
16 CFR 1633 and BFD IX-11

### Density and ILD

Top Layer – 3.0 lb. / 15 ILD

Top Layer Heel – 3.0 lb. / 15 ILD

Perimeter – 1.7 lb. / 55 ILD

Base Layer – 1.7 lb. / 30 ILD

*\*Mattress weight capacity only; total weight must not exceed bed frame manufacturer's specified load capacity.  
Joerns reserves the right to change specifications without notice.*

