



MANGAR Bathing Cushion

The Bathing Cushion is used to help someone who has difficulty getting in and out of the bath.

EASY TO FIT

- Fitted and removed within seconds
- Some air can be left in the backrest for greater levels of comfort
- The design adapts to the size and shape of the bath, so no expensive bathroom alterations are required

SIMPLE TO MAINTAIN

- Lightweight and portable, it can be removed and cleaned with ease
- Rinse with clean water after each bath and leave it to dry inflated
- Compact and easy to store
- The material is designed to exceptionally high standards with extraordinary durability, strength & performance in the toughest of environments

SAFE TO USE

- Only lowers a person into the bath if there is enough battery power to raise them back out again
- Battery operated – no electric cables
- Fully waterproof hand control (IP67)

IMPROVES WELL-BEING

- Perfect for a person who has difficulty getting in and out of the bath
- The cushion deflates allowing the user to lie flat in the bath, fully covered in water so benefiting from the therapeutic properties of bathing
- The user can take the bathing cushion away with them when travelling
- Can be removed from the bath when not in use so enabling a more flexible environment for shared bathrooms

TECHNICAL SPECIFICATIONS

Max user weight	330 lbs
Max height	15.75"
Min height	Flat
Seat depth	17.5"
Width	23"
Total product weight	4.5 lbs
Back/headrest length	13.75"
Airflo Compressor	12
Lifting performance per charge	Up to 8 lifts (4 bathing cycles)



Features

- 1. The air inlet valve can be rotated** allowing the hose to be led to the appropriate corner.
- 2. Pass the air hose through the loop in the corner of the backrest** on the side that you propose getting in and out of the bath.
- 3. Made with a durable, anti-microbial material.**
- 4. Strong suckers** ensure the bathing cushion doesn't move in the bath.
- 5. Safety feature** which ensures there is enough charge to lift user out of bath.
- 6. Seams joined by high frequency welding** for air tight construction




Accessories

- **Airflo 12 Compressor**



Product code: MPCA030520

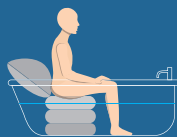

Joerns.com

General 800.826.0270

Canada 866.546.1151

6240084 RevA • 000576

US Partner with

