



MANGAR Surfer Bather

The Surfer Bather is a safe way of securely transferring a child in and out of the bath.

EASY TO FIT

- Fits into a bath by securing the base board suction cups firmly to the bath
- Elevates and lowers the child at the touch of a button
- Adjustable backrest and seat
- Adjustable pommel available

SIMPLE TO MAINTAIN

- Can be removed and cleaned with ease
- Rinse with clean water after each bath and leave it to dry inflated

SAFE TO USE

- Holds the child securely in the seat during the transfer into and out of the bath

- Includes side rails and leg straps for additional safety
- Fully waterproof hand control (IP67)
- Battery operated – no electric cables

IMPROVES WELL-BEING

- Perfect for parents and carers who struggle to lift a child into and out of the bath
- Bath time is a fun part of the day for both child and parent
- Lowers the child right to the bottom of the bath so they are completely covered in water and benefit from the therapeutic properties of bathing
- Remove from the bath when not in use to enable a more flexible environment for shared bathrooms

TECHNICAL SPECIFICATIONS

Max user weight	50 kg
Max user height	135 cm
Max height	48 cm
Min height	8.5 cm
Seat depth	62 cm
Total length	112 cm
Width (max-min)	67.5 cm 37 cm
Total product weight	8.7 kg
Backrest angle adjustment	Flat to 70°
Back/headrest length	56 cm
Airflo Compressor	12
Lifting performance per charge	Up to 20 lifts



Features


1. **Scissor lift** (including bellows and base board)
2. **Backrest** (including headrest pad)
3. **Platform pad with rails**
4. **Support straps**
5. **The backrest and the platform pad** are secured to the platform by means of plastic lugs inserted into the keyhold slots.
6. **Material is machine washable**

Accessories

- **Surfer Bather Platform Pad**
- **Pommel for Surfer Bather**
- **Head Support (Horseshoe)**
- **Surfer Bather Strap with Velcro**



Product code: MPCA040500

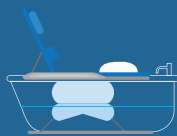

For more information visit

www.mangarhealth.com
 Email: enquiries@winn-care.uk
 Phone: 0800 2800 485

WIN  CARE